



## LEUKOCYTE ADSORBER LA 25


Adsorbent column for therapeutic leukocyte apheresis for IBD (inflammatory bowel disease) such as Ulcerative Colitis, Crohn Disease developed by MEDICA S.p.A, with cutting-edge features:

- Effective and Selective Removal of Activated Leukocytes and Platelets
- Excellent Biocompatibility thanks to **MediSulfone®** Membrane
- Excellent degree of safety and tolerability by patients
- Low Blood Volume
- Sterilization by Beta Radiation; No Gas Residual




## TEST IN VIVO

Blood parameters during each apheresis session  
(means of 10 patients; samples from leukocyte adsorber outflow)

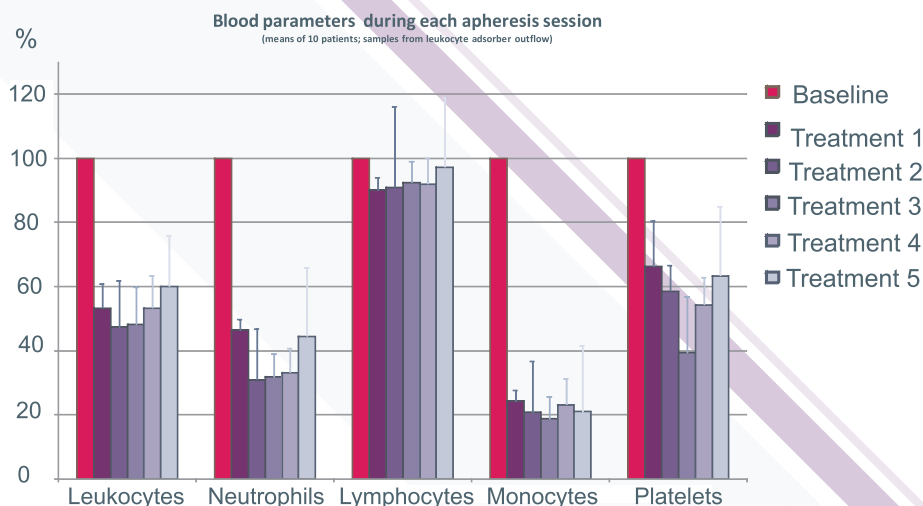

## SCHEMATIC DRAWING OF EXTRACORPOREAL CIRCULATION

LA 25 can be used on Medica AFERsmart Plus.

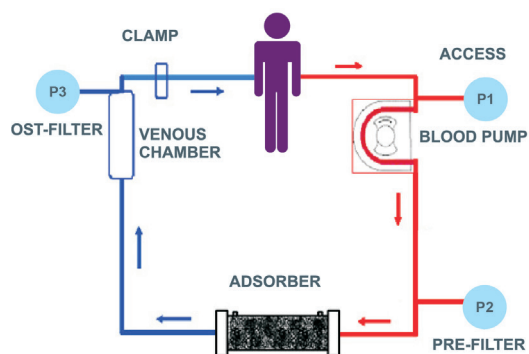

TECHNICAL DATA OF LEUKOCYTE ADSORBER LA25	
Fiber material	Polysulfone
Fiber length (mm)	MediSulfone®
Fiber surface (m <sup>2</sup> )	2
Fiber qty (gr)	25
Blood connectors	Twist Lock
Filter length (mm)	175
Treatment blood flow (ml/min)	30
Treatment duration (min)	60
Priming Volume	130 ml
Sterilisation	Beta Rays
Cartridges and blood port material	PC
Product code for Leukoadsorber kit	M03706

MediSulfone® is MEDICA's registered mark



Leukocyte adsorber LA 25 is assembled without latex, nor direct or indirect animal origin components.

### MANUFACTURER

MEDICA S.p.A. Via degli Artigiani, 7 | 41036, Medolla (Modena) Italia  
Tel: +39 0535 51159 | Fax: +39 0535 52605 | Email: info@medica-spa.com | www.medica-spa.com

**MEDICA**  
SEAMLESS INNOVATION FOR LIFE QUALITY

CE 0476 UNI EN ISO 9001:2015 UNI CEI EN ISO 13485:2021